


# Rapid Response™ Fentanyl (FYL) Test Strip



ONE LINE  
FENTANYL

TWO LINES  
NO FENTANYL

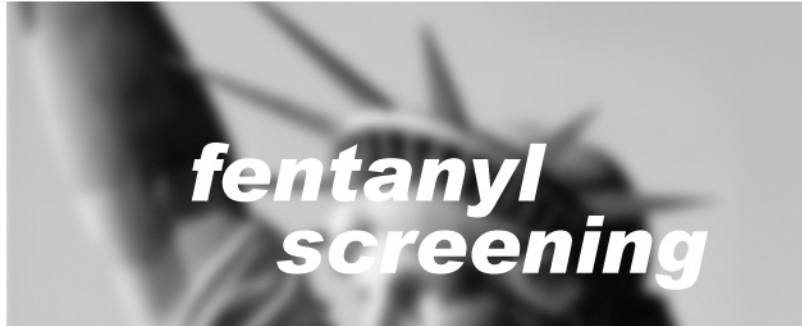

The Fentanyl Forensic Test Kit is a lateral flow Chromatographic Immunoassay for the qualitative detection of Fentanyl in liquid and powder substances at the cut-off concentration of 200 ng/mL.

Both are specific for Fentanyl screening with no significant cross reactivity to other opiates, such as Morphine and Heroin. The Rapid Response™ Fentanyl (FYL) Forensic Test Kit is able to detect Fentanyl and many other Fentanyl analogues such as Carfentanil, Acetyl Fentanyl, Butyryl Fentanyl, Remifentanyl, Ocfentanil, Sufentanil, p-Fluoro Fentanyl, Furanyl Fentanyl, Valeryl Fentanyl, and 3-Methyl Fentanyl.

Interpretation of test results are: positive (one line), negative (two lines), invalid (no lines or no control line). Test results are available in 5 minutes after test completion.

## Forensic Test Strip/s

### Product Information – Fentanyl 100 Kit

Product Code:	LS_Fentanyl100
Item #:	FYL-18S7-100
Sample:	Liquid / Powder
Format:	Strip
Quantity:	100 Tests/Kit
Sensitivity	200 ng/mL
Specificity:	Fentanyl
Accuracy:	> 98%
Time-to-Result:	5 minutes
Storage Condition:	2-30°C/36-86°F
Test Principle:	Lateral Flow Immunoassay
Contents:	100 Test strips, Package Insert, 100 Vials Proprietary Buffer Solution (2 ml/ea)

### Product Information – 5 Test Kit

Product Code:	LS_FentanylIRRK
Item #:	FYL-18S7-5
Sample:	Liquid / Powder
Format:	Strip
Quantity:	5 Tests/Kit
Sensitivity	200 ng/mL
Specificity:	Fentanyl
Accuracy:	> 98%
Time-to-Result:	5 minutes
Storage Condition:	2-30°C/36-86°F
Test Principle:	Lateral Flow Immunoassay
Contents:	Test strips, Swabs, Test Tubes, Extraction Buffers, Workstation,

